
8 + DHS

80 + DHS

pH - Cond - PC

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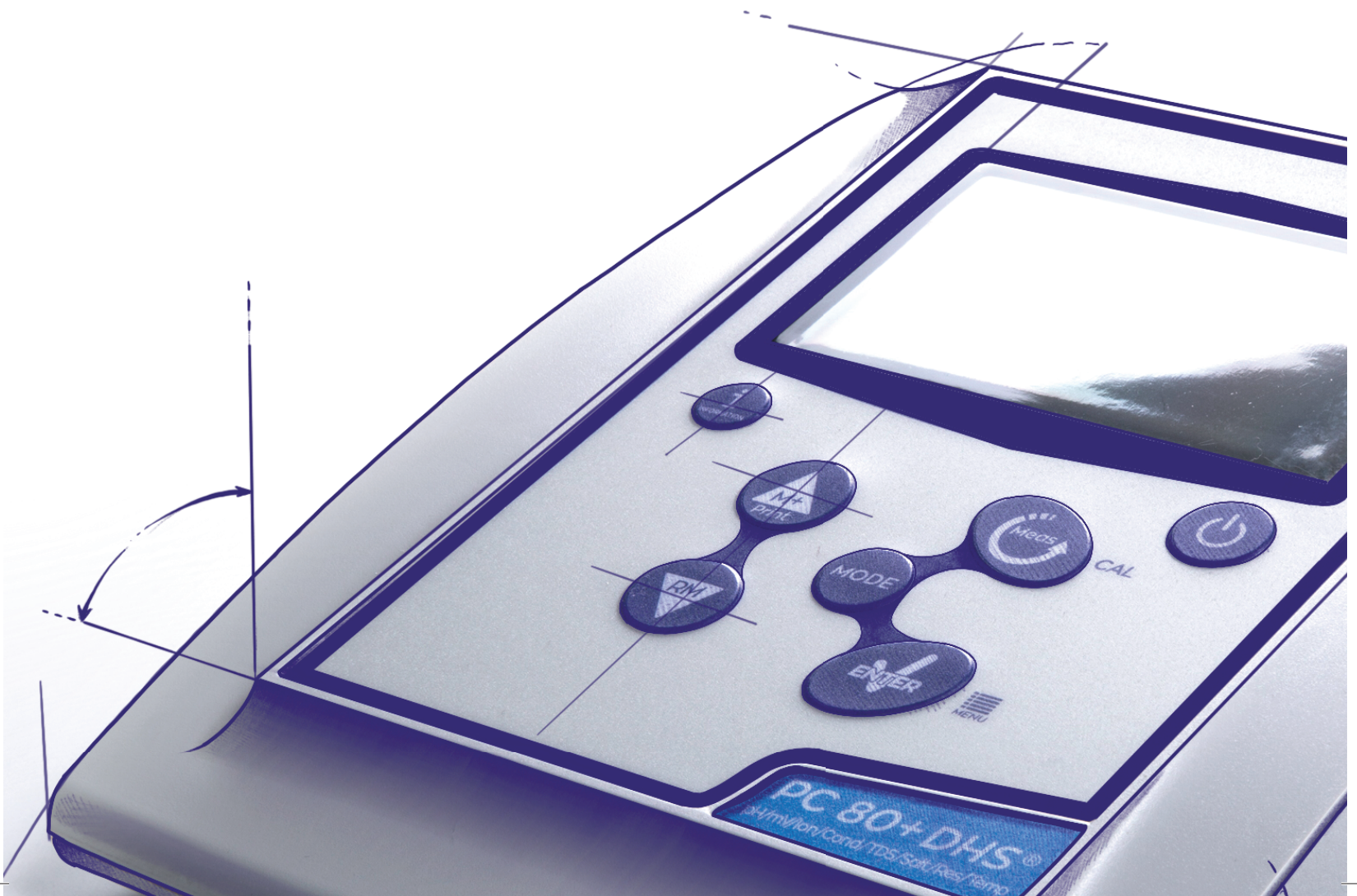INDEX

1. Introduction	4
2. Safety information	5
• Definition of warning words and symbols	5
• Reporting terms	5
• Additional documents for safety	6
• Use according to destination	6
• Basic requirements for a safe use	6
• Unauthorized use	6
• Device Maintenance	6
• Responsibility of the owner of the instrument	7
3. Instrumental Features	7
• Parameters	7
• Datasheet	7
4. Description of the instruments	8
• Display	8
• Keyboard	9
5. Installation	9
• Supplied components	9
• Start up	9
• Connection of the power supply	9
• Turning on and off, update date, time and language	10
• Instrument transportation	10
• Key functions	11
• Inputs / Outputs connections	11
• Symbols and icons on the display	12
6. Operation of the device	12
7. Setup Menu	13
8. Temperature Measurement ATC – MTC	13
9. pH parameters	14
• pH parameters setup	14
• Automatic pH calibration	16
• Calibration with manual values	17
• Performing pH measurement	18
• DHS Electrode	18
• Errors during calibration	19
10. ORP Parameter (Redox Potential)	19
11. Measurement with Ion-Selective Electrodes (ISE/ION)	19

• Ion parameters setup	19
• Calibration with Ion-Selective Electrodes	21
• Measurement with Ion-Selective Electrodes	22
12. Conductivity Parameter	22
• ... how to get to the conductivity?	22
• Configuration of setup menu for Conductivity parameter	23
• Automatic Conductivity calibration	25
• Manual conductivity calibration	26
• Performing a Conductivity measurement	27
13. Other measurements perform with conductivity cell	28
• TDS	28
• Salinity	29
• Resistivity	29
14. Data Logger and printer functions	29
• Setup for Data Logger parameter	29
• Example 1 Data Logger mode	30
• Example 2 Data Logger mode	30
• Recalling saved data on instrumental memory	31
• How to delete data in memory	31
15. Configuration menu	31
16. Magnetic stirrer with independent control	33
• Functioning	33
17. Software DataLink+ (for Windows 7/8/XP/10)	33
• Functions	34
18. Warranty	35
19. Disposal of electronic devices	35

1. Introduction

XS Instruments, globally recognized as a leading brand in the field of electrochemical measurements, has developed this new line of bench instruments completely produced in Italy, finding the perfect balance between performance, attractive design and ease of use. The large, simple and intuitive color display shows all the necessary information such as the measurement, the temperature, the buffers used for the last calibration (also custom), the electrode conditions, GLP information and stored data.

Everyone can use these tools thanks to the instructions that appear directly on the display. The interface is multilingual (8 different languages available) and the English operation guide for Series 80+ is available at the push of a button.

Up to 5 pH calibration points (3 for Series 8+) can be performed using the USA NIST and DIN families with automatic recognition and 4 points for conductivity; also buffers chosen by the operator can be used.

A calibration can be saved for the conductivity measurement for each cell constant used; the compensation coefficient for ultra pure water is also available in the 80+ series.

A detailed calibration report with innovative graphic representation and the representation through the icons of the buffers used make the calibration process more efficient.

For the 80+ series, the measurement of Selective Ions is available with standards and units of measurement that can be set by the user. Automatic or manual Data Logger with values that can be stored in different GLP formats on the internal memory (1000 data), on the PC or printable in paper format.

Password for calibration management, memory emptying and date / time modification. The ideal solution for an accurate and precise measurement is to use an XS Instruments electrochemical sensor from the XS Sensor range with an XS Instruments device and perform the calibrations by supplying XS Solution certified calibration solutions.

2.Safety information

- **Definition of warning words and symbols**

This manual contains extremely important safety information, in order to avoid personal injury, damage to the instrument, malfunctions or incorrect results due to failure to comply with them. Read entirely and carefully this manual and be sure to familiarize with the tool before starting to work with it.

This manual must be kept near to the instrument, so that the operator can consult it easily, if necessary. Safety provisions are indicated with warning terms or symbols.

- **Reporting terms:**

ATTENTION for a medium-risk hazardous situation, which could lead to serious injury or death, if not avoided.

ATTENTION for a dangerous situation with reduced risk which can cause material damage, data loss or minor or medium-sized accidents, if not avoided.

WARNING for important information about the product.

NOTE for useful information about the product.

Warning Symbols:

Attention

This symbol indicates a potential risk and warns you to proceed with caution.

Attention

This symbol draws attention to a possible danger **from electric current**.

Attention

The instrument must be used following the indications of the reference manual. Read the instructions carefully.

Advice

This symbol draws attention to possible damage to the instrument or instrumental parts.

Note

This symbol highlights further information and tips.

• Additional documents for safety

Following documents provide the operator with additional information to work safely with the measurement system:

- Operating manual for electrochemical sensors
- Safety data sheets for buffer solutions and other maintenance solutions (e.g. storage)
- Specific notes on product safety.

• Use according to destination

This instrument is designed exclusively for electrochemical measurements both in laboratory. Pay attention to the technical specifications shown in the INSTRUMENT FEATURES / TECHNICAL DATA table; any other use is to be considered unauthorized. This instrument has left the factory in perfect technical and safety conditions (see test report in each package). The regular functionality of the device and the operator safety are guaranteed only if all the normal laboratory safety standards are respected and if all the specific safety measures described in this manual are observed. All the notes, indications and recommendations contained in this manual are also valid for the independently controlled magnetic stirrer that can be found in combination with some sales codes for the instruments.

• Basic requirements for a safe use

The regular functionality of the device and the operator safety are guaranteed only if all the following indications are respected.

- The instrument can be used in accordance with the specifications mentioned above only.
- Use the supplied power supply only. If you need to replace the power supply, contact your local distributor.
- The instrument must operate exclusively in the environmental conditions indicated in this manual.
- Neither the instrument nor the magnetic stirrer should be opened by the user for any reason.

Do this only if explicitly authorized by the manufacturer.

• Unauthorized use

The instrument should not be put into operation if:

- It is clearly damaged (for example due to transportation);
- It has been stored for a long period of time in adverse conditions (exposure to direct light, heat sources or places saturated by gas or vapours) or in environments with conditions different from those mentioned in this manual.

• Device Maintenance

If used correctly and in a suitable environment, the instrument does not require maintenance procedures. It is recommended to occasionally clean the instrument case with a damp cloth and a mild detergent. This operation must be performed with the instrument off, disconnected from the power supply and by authorized personnel only. The housing is in ABS / PC (acrylonitrile butadiene styrene / polycarbonate). This material is sensitive to some organic solvents, for example toluene, xylene and methyl ethyl ketone (MEK). If liquids get into the housing, they could damage the instrument. In case of prolonged non-use of the device, cover the BNC connectors with the special cap supplied with. Do not open the instrument housing: it does not contain parts that can be maintained, repaired or replaced by the user. In case of problems with the instrument, contact your local distributor. It is recommended to use original spare parts only. Contact your local distributor for information. The use of non-original spare parts can lead to malfunction or permanent damage to the instrument. Moreover, the use of spare parts not guaranteed by the supplier can be dangerous for the user himself. For the maintenance of the electrochemical sensors, refer to the documentation present in their packaging or contact the supplier.

• **Responsibility of the owner of the instrument**

The person who owns and uses the device or authorizes its use by other people is the owner of the instrument and is responsible for the safety of all users of the device and third parties

The owner of the instrument must inform users of the use of the same safely in their workplace and on the management of potential risks, also providing the required protective devices

When using chemicals or solvents, follow the manufacturer's safety data sheets.

3. Instrumental Features

• **Parameters**

pH8+DHS: pH, ORP, Temp

Cond8+: Cond, TDS, Sal, Temp

PC8+DHS: pH, ORP, Cond, TDS, Sal, Temp

pH80+DHS: pH, ORP, Ion, Temp

Cond80+: Cond, TDS, Sal, Res, Temp

PC80+DHS: pH, ORP, Ion, Cond, TDS, Sal, Res, Temp

• **Datasheet**


	8+ Series	80+ Series
pH	pH8+DHS - PC8+DHS	pH80+DHS - PC80+DHS
Measuring Range	-2...16	-2...20
Resolution / Accuracy	0.1, 0.01 / ±0.02	0.1, 0.01, 0.001 / ±0.002
Calibration points and buffers recognition	1...3 / USA, NIST, 2 user values	1...5 / USA, NIST, DIN, 5 user values
Indication of calibration points	Yes	
Recognition DHS electrode	Yes	
pH alarm threshold MIN/MAX	Yes	
Analogic display	Yes	
Stability filter	Med - High - Tit	
ORP	pH8+DHS - PC8+DHS	pH80+DHS - PC80+DHS
Range / Resolution	±2000 / 1 mV	±2000 / 0.1, 1 mV
ISE	pH8+DHS - PC8+DHS	pH80+DHS - PC80+DHS
Resolution	-	0.001 - 0.099 / 0.1 - 19.9 / 20 - 199 / 200 - 19999
Calibration points	-	2...5
Unit of measure	-	mg/L - g/L - mol/L
Conductivity	Cond8+ - PC8+DHS	Cond80+ - PC80+DHS
Range / Resolution	0,00 - 20,00 - 200,0 - 2000 µS / 2,00 - 20,00 - 200,0 mS Automatic scale	0,00 - 20,00 - 200,0 - 2000 µS / 2,00 - 20,00 - 200,0 - 1000 mS Automatic scale
Calibration points and buffers recognition	1...4 / 84, 1413 µS, 12.88, 111.8 mS, 1 user values	
Refernce temperature	15...30 °C	
Temperature compensation factor	0,00...10,00 %/°C	0,00...10,00 %/°C and Ultrapure Water
TDS	Cond8+ - PC8+DHS	Cond80+ - PC80+DHS
Measuring Range / TDS Fatctor	0,1mg/L...100 gr/L / 0.40...1.00	0,1mg/L...500 gr/L / 0.40...1.00
Salinity	Cond8+ - PC8+DHS	Cond80+ - PC80+DHS
Measuring Range	0,01...100ppt	

Resistivity	Cond8+ - PC8+DHS	Cond80+ - PC80+DHS
Measuring Range	-	0....10 MΩ*cm
Temperature		
Measuring Range	-10...110 °C	-20...120 °C
Resolution / Accuracy	0,1 / ± 0,2°C	
Temperature compensation ATC (NTC30KΩ) and MTC	0...100 °C	
System		
GLP with Calibration Timer	Yes	
Memory	1000 Data	
Display	Multicolor dot matrix	
Multilingual	Yes	
Quick On-line manual	-	Yes
IP protection	IP 54	
Power Supply	Power Supply AC/DC 12V / 1000mA	
Tolerance relating to the power supply	± 10%	
Electrical voltage	100 – 240 V	
System		
Work frequency	47 – 63 Hz	
Maximum absorption	10 W	
Sound level during standard operation	< 80 dB	
Environmental storage conditions	-25 ... +65 °C	
Environmental operating conditions	0 ... +45 °C	
Maximum permissible humidity	< 95 % non-condensing	
Grade of micro-pollution of the environment in which the product is used	grade 2	
Maximum altitude of use	2000 m	
System dimensions LxPxA	Vers. Basic: 162 x 185 x 56 mm	Vers. Stirrer: 324 x 185 x 56 mm
System weight	Vers. Basic: 709 g	Vers. Stirrer: 1255 g


4. Description of the instruments

• Display

* function available only for 80+ Series

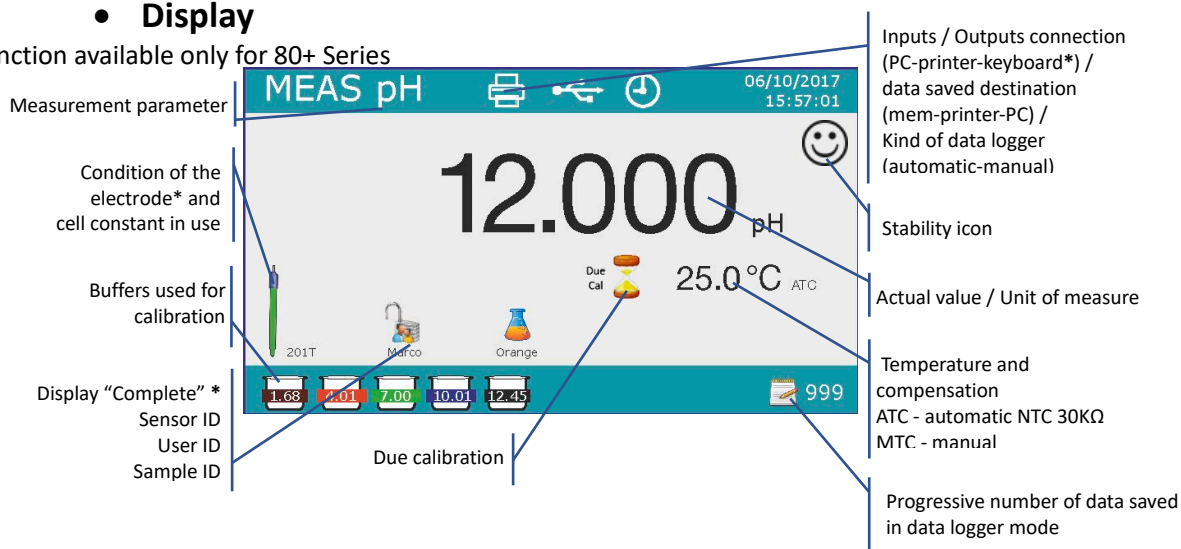• Keyboard

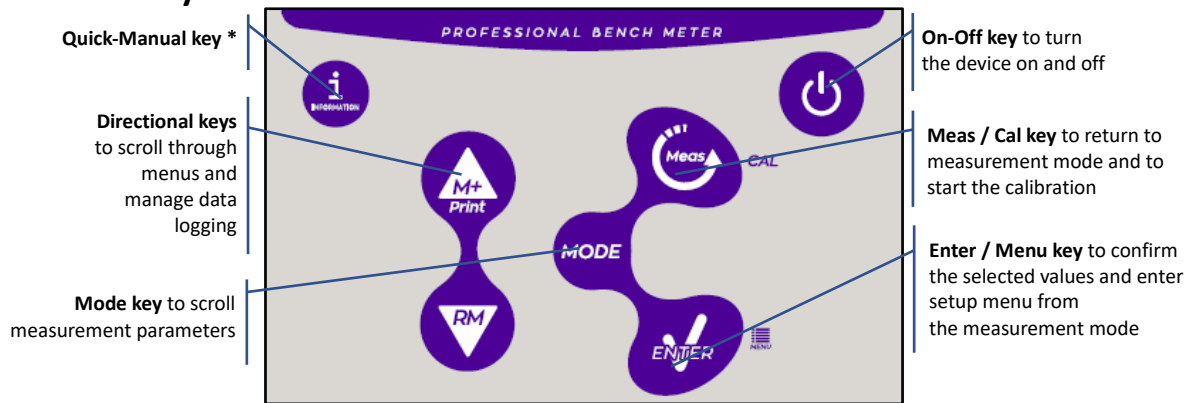8+ Series keys

80+ Series keys

5. Installation

• Supplied components

In the BASIC versions, it is always included inside the package:

Instrument, 12V / 1000mA power supply, electrode holder, 1m S7 / BNC connection cable, NT55 temperature probe, single-dose bottle buffer solutions, USB cable, (external keyboard only for 80+ series), multilingual user manual and test report.

In the STIRRER versions, it is always included inside the package:

Instrument, 12V / 1000mA power supply, independently controlled magnetic stirrer with connection cable, metal anchors, electrode holder arm, 1m S7 / BNC connection cable, NT55 temperature probe, buffer solutions in single-dose bottle, USB cable, (external keyboard only for 80+ series), multilingual user manual and test report

• Start up

- Place the instrument on a flat, stable laboratory bench with adequate front and side accessibility. It is recommended to place the instrument at a distance of not less than 20 cm from overlying and surrounding parts.
The positioning carried out in this way eliminates the residual risk of possible minor damages caused by manual handling of the loads.
- Make sure that the instrument and the surrounding space are correctly illuminated
- In the STIRRER versions, using a Phillips screwdriver, unscrew the metal brackets located under the stirrer and attach them to the device. Then insert the electrode holder arm on the pin. Connect the instrument ("Stirrer" connector) at the stirrer with the special cable provided.
- In the BASIC versions, assemble the holder next to the instrument by assembling the electrode holder arm into the base pin.
- The instrument is not designed to be used in atmospheres where hazardous materials are present.

• Connection of the power supply

- Check that the electrical standards of the line on which the instrumentation is to be installed comply with the voltage and operating frequency of the instrument (technical datasheet).**
- Use the original power supply only.**
- Connect the power supply plug to the connector on the rear panel of the instrument indicated with the icon $12\text{ V} \text{ --- } \text{---}$
the icon
- Connect the power supply to an electric socket easy to reach.

- The instrument is equipped with an external power supply that has no protection against infiltration of liquids, therefore for its use it is necessary to keep all the electrical cables and connections away from any liquids and humidity and not to use the appliance in a humid room such as a bathroom or laundry room.

ATTENTION - Danger of death or serious injury from electric shock. 

Contact with live components can lead to injury or death.

- Use only the adapter supplied.
- Do not put the power supply in contact with liquids nor in a condensing environment, avoid thermal shock.
- All electrical cables and connections must be kept away from moisture or liquids.
- Check that all cables and plugs are not damaged, otherwise replace them.
- During use, do not cover the power supply and/or do not place it inside boxes
- In the event of a loss of power while the appliance is operating, there is no dangerous condition when the appliance is powered again. In this case, it will be necessary to restart it as it will not restart automatically.

• Turning on and off, update date, time and language

Connect the instrument to the electricity grid using only the supplied power supply, press  key, the following screens appear on the display:

- Device model and software.
- Settings for the most important parameters and any info on DHS sensor.

At first use it is advisable to update the date, time and language of the device by performing the following operations:

- Press  key to enter in menu setup.
- Press  up to select " Configuration" and confirm with .
- Press  up to select **P6.10 Date Setting** and confirm with .
- Use  and  keys to update the "year" value, press  to confirm and move to "month", repeat this operation for "day" too; the device automatically updates the date and returns to the setup menu.
- Press  , select **P6.11 Time Setting** and confirm with  , use  and  keys and confirm with  to update "hour", "minutes" and "seconds".
- Press twice  key, select **P6.13 Select Language** and confirm with  , move with  on the desired language and activate it with .
- To turn off the device press  key in measurement mode.

• Instrument transportation

To move the instrument to a new location, transport it carefully to avoid damage; the instrument can be damaged if it is not transported correctly. Disconnect the instrument from the power supply and remove all the connection cables. Remove the electrode arm. To avoid damage to the instrument during long distance transport, use the original packaging. If the original packaging is no longer available, choose a package that guarantees safe transport.

• **Key functions**

Key	Push	Functions
	Short	Press to turn the device on or off.
	Short	In measurement mode press to scroll through the parameters: <ul style="list-style-type: none"> • pH80+DHS: pH → ORP → analogic pH → Ion* • Cond80+: Cond → TDS → Sal → Res* • PC80+DHS: pH/Cond → pH → ORP → analogic pH → Ion* → Cond → TDS → Sal → Res* * parameters available only for 80+ Series
	Short	In calibration, setup and recall memory mode press to return in measurement mode.
	Long (3s)	In measurement mode press and hold to start the calibration.
	Short	In measurement mode press to enter in setup menu. In setup menus, press to select the program and / or the desired value. During calibration, press to confirm the value.
	Short	In setup and subsetup menu press to scroll. In the submenu of the setup press to change the value. In memory recall mode press to scroll the saved values. In MTC mode and customer calibration press to change the value. M+/Print: In measurement mode press to save or print the value (maual data logger) or start and finish the recording (automatic data logger). RM: In measurement mode press to recall saved data.
		Long (3s)
	Short	In measurement press to display the quick manual with the calibration instructions.


The correct use of the function keys and the attention in pressing them, given the small size, eliminates the residual risk of minor damages, not probable, caused by simultaneous pressing of the keys; before each use, check that pressing the keys has the corresponding effect on the display.

• **Inputs / Outputs connections**

Use only original accessories guaranteed by the manufacturer.

If necessary, contact your local distributor. The BNC connectors at the time of sale are protected by a plastic cap. Remove the cap before connecting the probes.

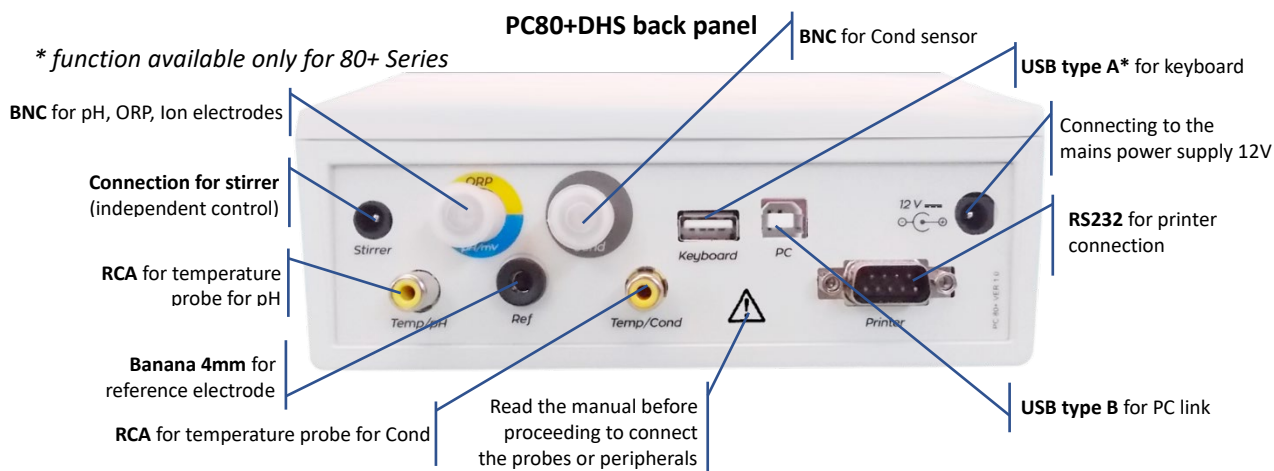• **Symbols and icons on the display**


Symbol	Description	Symbol	Description
	Number of data stored in Data Logger mode on instrument memory		Titration stability criterion (continuous measurement)
	Printer data logger		Measurement stability icon
	External keyboard*	HOLD	HOLD mode, value is frozen
	Data logger with data sending on PC		No user Password
	Data logger with data sending on instrument memory		User password entered
	Automatic data logger mode (when blinking is recording)		Sample ID*
	Manual data logger mode		Calibration expired

* function available only for 80+ Series


6. Operation of the device

- When the device is powered on, the instrument enters measurement mode in the last parameter that was used.
- To scroll through the parameter screens, press key, the current measurement parameter is shown on the upper left display.

pH8+DHS - pH80+DHS


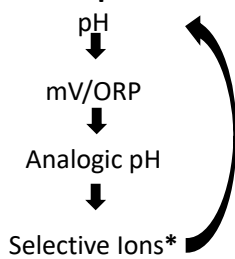
COND8+ - COND80+

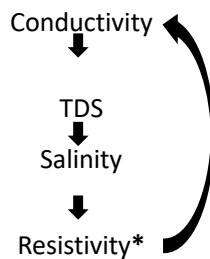PC8+DHS - PC80+DHS

Simultaneous pH/Cond

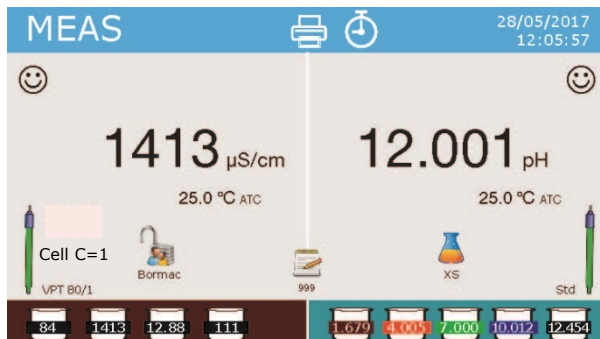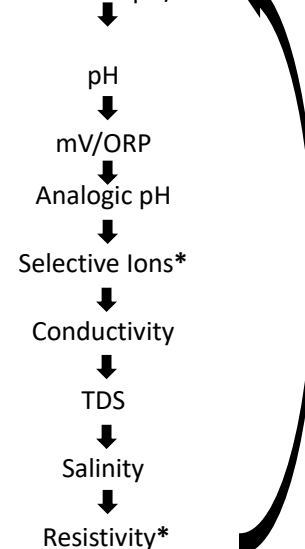Example: simultaneous screen pH/Cond on PC8+DHS and PC80+DHS

* parameters available only for 80+ Series

- In the measurement screen you are using press and hold "CAL" key to start calibration of the active parameter (see next paragraphs).


Reading this manual before each use eliminates the residual risk of possible and significant errors in the interpretation of the data on the display. It is therefore recommended that the user carefully read the manual to correctly use the instrument and correctly interpret the information on the display. All this will ensure that the risk of misinterpretation goes from possible to highly improbable.

7.Setup Menu

- In measurement mode press “ MENU” key to enter in SETUP mode, choose the parameter that you want to change by moving with the directional keys and confirm with .
- pH8+DHS:** pH ↔ Data Logger ↔ Configuration
- pH80+DHS:** pH ↔ Ion ↔ Data Logger ↔ Configuration
- Cond8+ and Cond80+:** Cond/TDS ↔ Data Logger ↔ Configuration
- PC8+DHS:** pH ↔ Cond/TDS ↔ Data Logger ↔ Configuration
- PC80+DHS:** pH ↔ Ion ↔ Cond/TDS ↔ Data Logger ↔ Configuration
- In the selected menu, move between different programs using directional keys and press key to select what you want to change.
- Using and keys choose the option you want or change the numeric value and confirm with key.
- Press key to return in measurement mode.

Note: External keyboard for 80+ Series can be used to scroll through the different setup menus

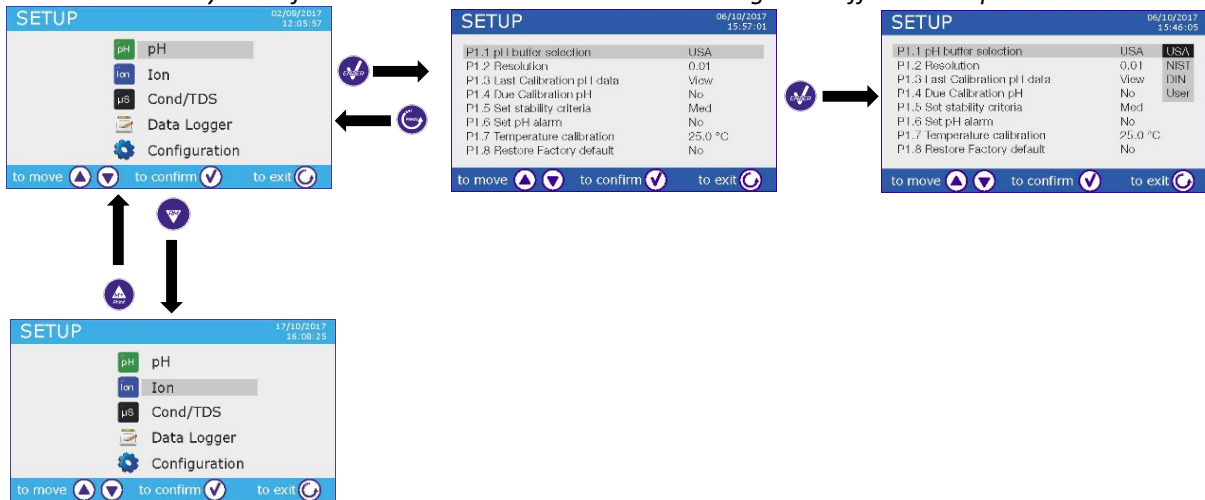8.Temperature Measurement ATC – MTC

- ATC:** The direct measurement of the sample temperature for all parameters is carried out through the NTC 30KΩ probe, which can either be integrated into the electrode or external.
- MTC:** If no temperature probe is connected, this must be manually inserted: press and hold or until two little arrows appear alongside the value, adjust the value using again directional keys, then press key to confirm.

9.pH parameters

pH8+DHS, pH80+DHS, PC8+DHS, PC80+DHS

On this series of devices, you can use pH sensors with integrated temperature probe or connect two different sensors. The pH electrode uses a BNC connector while the temperature probe uses a RCA/CINCH connector. The instrument is also able to recognize DHS sensor, an innovative electrode capable of storing the calibration data and then being immediately used on any enabled instrument.

• pH parameters setup

- In measurement mode press “ MENU” key to enter SETUP menu.
- Use  and  keys to scroll until “pH” menu and press  key to confirm.
- Use  and  keys to move and then select the program you want to edit.

The table below shows the structure of the setup menu for the pH parameter, for each program you can find the options that the user can choose and the default value.

Program	Description	Options	Factory default
P1.1	pH buffer selection	USA – NIST – DIN* – USER	USA
P1.2	Resolution pH	0.1 – 0.01 – 0.001*	0.01
P1.3	Last calibration pH data	View data – Print	View data
P1.4	Due calibration pH	NO – 1...99h – 1...99d	NO
P1.5	Set stability criteria	Tit – Med - High	Med
P1.6	Set pH alarm	NO – MIN - MAX	NO
P1.7	Temperature calibration	YES – NO	-
P1.8	Restore factory default	YES – NO	NO

* function available only for 80+ Series

P1.1 pH buffer selection

- Access this setup menu to choose the family buffer for calibrating the electrode.
- The **80+ Series** device allows the pH calibration to be performed from **1 to 5 points**, while **8+Series** can perform calibrations from **1 to 3 points**.
- During the calibration press  key to exit and save the calibrated points up to that moment (see calibration chart).
- The instrument automatically recognizes 3 buffers family (**USA, NIST** and **DIN***) in addition, the user has the possibility to calibrate manually up to 5 points with **customizable** values (2 points for 8+ Series).

USA buffers: 1,68 - 4,01 - **7,00**** - 10,01 - 12,45 (factory default)

NIST buffers: 1,68 - 4,00 - **6,86**** - 9,18 - 12,46

DIN buffers: 1,68 - 4,01 - **6,86**** - 9,18 - 12,45

** The neutral point is always required as the first point

In measurement mode at the bottom left of display a set of becker indicates the tampons with which the last calibration was made both automatically and manually. Inside the becker the number represents the exact value of the buffer, and for a quick and intuitive understanding a color scale has been inserted.

Becker color	Buffer pH value
Brown	< 2.5
Red	2.5 ~ 6.5
Green	6.5 ~ 7.5
Blue	7.5 ~ 11.5
Black	> 11.5

P1.2 Resolution pH


Access this menu to choose the resolution you want to get in measurement and calibration pH of the parameter:

- **0.1**
- **0.01** -factory default-
- **0.001** (available only for **80+ Series**)

P1.3 Last calibration pH data

Access this menu to get information on the last performed calibration.

- “View” (Factory default), a report is shown on the display with the following information about the current calibration:
CALIBRATION DATE and HOUR / TEMPERATURE / DHS DATA / OFFSET / SLOPE % for every range.
- Press  key to get to the innovative **graphical representation** of the calibration conditions so you can intuitively view the status of the sensor.

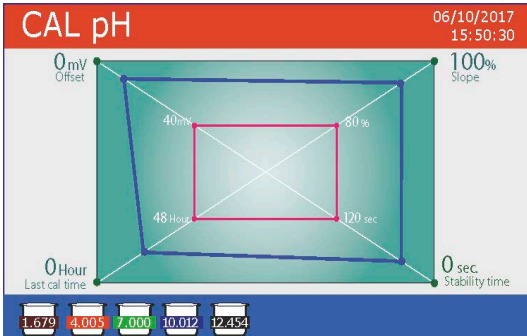The graphical calibration report is designed to provide the user with an immediate view of the calibration conditions, indeed, the closer are the blue lines (calibration data) to the outside of the graph, the closer you are to the ideality of calibration and electrode conditions; conversely, the condition worsens the more you approach the red rectangle which represents the limit of acceptability. The chart shows offset data, average slope and sensor stabilization time, and how many hours have passed since the last calibration.

- “Print”: connect the printer to the device through RS232 port (see outputs connections) to print the calibration report directly on paper format.

P1.4 Due calibration pH

Enter this menu to set a calibration expiration, this option is essential in GLP protocol.

- No calibration expiration is set by default, use  and  keys to choose days or hours which must pass between two calibrations and confirm with  key.
- Once the calibration has expired, “” icon will appear on the display, anyway the user may continue to make measurements.

Note: The parameter “Last calibration time” in the graphical calibration report is not affected by the expired calibration, it is just informative.

Note2: With DHS electrode the expired calibration is related to the electrode.

P1.5 Set stability criteria

To consider the reading of a value truthful it is advisable to wait for measurement stability, which is indicated by the icon .

Enter in this setup menu to set stability criteria.

- “Medium” (factory default): measurements included within 0.6 mV.
- “High”: choose this option for a more rigorous reading, measurements included within 0.3 mV.
- “Tit” ((titration) no stability criteria is activated, the measurement will be “continuous”.

With this option activated, the icon  appears on the display and the measurement will hardly stabilize, however, the instrument response time is minimized as it is a simultaneous measurement.

P1.6 Set pH alarm

The user can set an alarm in the **Low and / or Maximum pH value**.

When the threshold value is exceeded, the device will emit a beep. In the analog display mode, the alarm range is indicated by the red color.

P1.7 Temperature calibration

All the instruments within this series are pre-calibrated for a correct measure of the temperature. However, if there is a difference between the measured and the real one (usually due to a malfunction of the probe), it is possible to adjust the offset of + 5 ° C.

- Use  and  keys to correct the temperature offset value and press  key to confirm.

P1.8 Restore factory default.

If the instrument does not properly work or incorrect calibration were performed press **Yes** with  key to return all pH parameters to the factory default settings.

IMPORTANT: the factory reset of the parameters does not delete stored data.

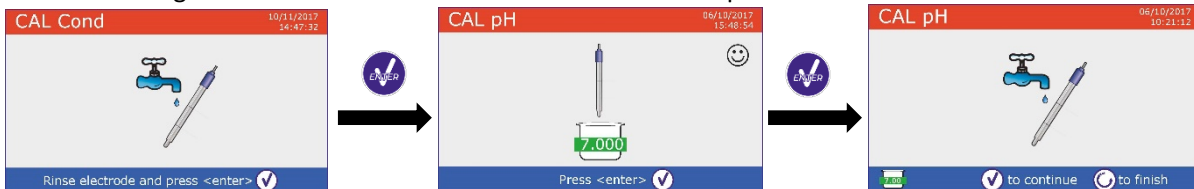
• Automatic pH calibration

Example for 3 points USA buffers calibration


- In pH measurement mode press and hold for 3 sec “ CAL” key to enter into calibration mode.
In PC devices pH calibration can be accessed also from simultaneous pH/Cond screen then selecting **pH**.

- Rinse the electrode with distilled water  and gently dab with absorbent paper.
Press  and immerse the electrode in the pH buffer 7.00 (as indicated by the becker on the display). The first calibration point is always the neutral pH (7.00 for USA buffer, 6.86 for NIST and DIN buffers) while the others are at the discretion of the user.

- When the  icon appears on the display confirm the first calibration point pressing  key. The display flashes the value measured and then the pH 7.00 icon becker appears at the bottom left  indicating that the instrument is calibrated at the neutral point.

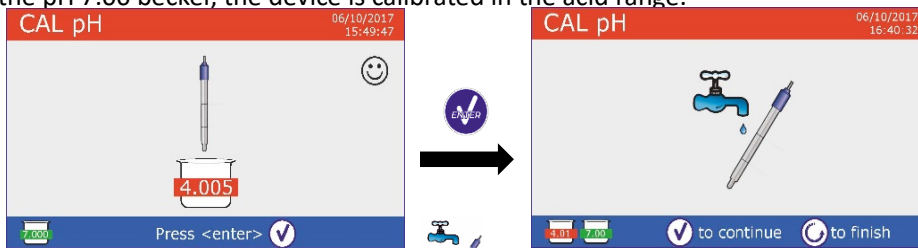- Extract the electrode, rinse with distilled water  and gently dab with absorbent paper.
- Press  key to continue the calibration and immerse the electrode in the pH buffer 4.01. In the becker all the values that the device can recognize flow.
- When the value on the becker stabilizes on pH 4.01 and the  icon appears confirm by pressing  key. The display flashes the value measured and subsequently the pH 4.01  icon becker appears alongside the pH 7.00 becker, the device is calibrated in the acid range.

- Extract the electrode, rinse with distilled water  and gently dab with absorbent paper.
- Press  key to continue the calibration and immerse the electrode in the pH buffer 10.01. In the becker all the values that the device is able to recognize flow.
- When value on the becker stabilizes on pH10.01 and the  icon appears confirm by pressing  key.


The passage from an acid to a basic pH may take a few more seconds to achieve stability.

- The display flashes the value measured and subsequently the pH 10.01 icon becker  appears alongside the pH 7.00 and pH 4.01 icon becker, the device is calibrated also in alkaline range.
- Although the device could accept two other calibration points, we stop and confirm this calibration pressing  (**8+ Series** will automatically end after three points).

The calibration report and graphic representation are displayed, press  or  to exit and return in measurement mode. Bottom left are displayed buffers used for the last calibration.

- Example for 3 Points Calibration for 8+ Series

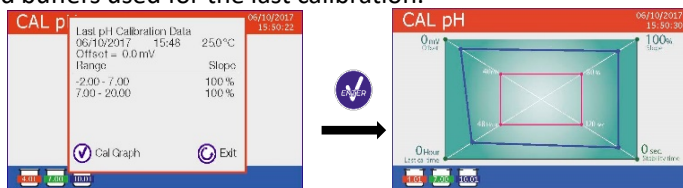- Example for 5 Points Calibration for 80+ Series

Note: calibration of the electrode is a fundamental operation for the quality and veracity of measurement. Make sure that the buffers used are new, not polluted and at the same temperature. After a long time or after reading special samples renew the calibration, the graphical report can help the user making this decision.

ATTENTION: Before proceeding with the calibration operations, carefully consult the safety data sheets of the substances involved:

- Calibration buffer solutions.
- Storage solution for pH electrodes.
- Filling solution for pH electrodes.

The careful reading of the safety data sheets of the solutions used favors the elimination of residual risks related to skin, contact, ingestion, inhalation or eye contact that can generate possible but not probable minor damages. Contact your local distributor for more information.

• Calibration with manual values

Example for a two-points calibration pH 6.79 and pH 4.65 (DIN19267)

- Enter in the pH setup menu and set in P1.1 → USER, press  to return in measurement pH mode.
- Press and hold for 3 sec “ CAL” key to enter in calibration mode.


- Rinse the electrode with distilled water  and gently dab with absorbent paper.

- Press  and immerse the electrode in the first buffer solution (ex pH 6.79).

- Wait for the stability, when  icon appears use  and  keys to adjust the value by entering the real buffer value (ex pH 6.79).

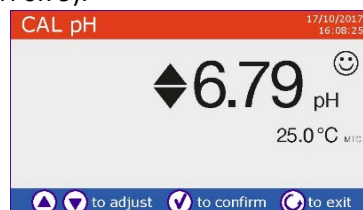
Note: Check the buffer value according to the temperature

- When  icon reappears press  key to confirm the first point, the actual measured value flashes on the display and the becker icon appears with the identifier color and the buffer value at the bottom left  U 6.79 (U means “USER value”).

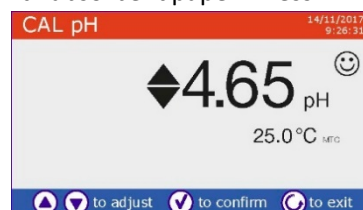- Extract the electrode, rinse with distilled water  and gently dab with absorbent paper. Press  key to continue calibration and immerse the electrode in the next buffer (ex pH 4.65).

- Wait for the stability, when  icon appears use  and  keys to adjust the value by entering the real buffer value (es pH 4.65).

- When ☺ icon reappears press  key to confirm the second point, the actual measured value flashes on the display and the becker icon appears with the identifier color and the buffer value at the bottom left .
- Although the device could accept two other calibration points, we stop and confirm this calibration pressing  (**8+ Series** will automatically end after three points).
- The calibration report and graphic representation are displayed, press  or  to exit and return in measurement mode. At bottom left the buffers used for the last calibration are displayed, the value is preceded by the letter "U" indicating that the value was entered manually. 


Note: If you are working with manual temperature compensation (MTC), update the value before calibrating the device.

• Performing pH measurement


- Enter in the Setup menu for pH to check the calibration and if necessary, update the reading parameters (see paragraph “pH parameters setup”), press  to return in measurement mode.
- Press  key to scroll through the different parameter screens until you find **MEAS pH** (see paragraph “Operation of the device”).

The pH measurement can be displayed on three different screens:


• **Digital**

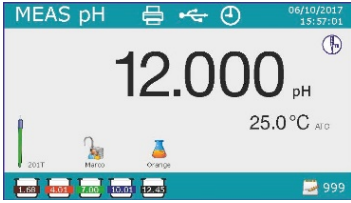• **Analogic**

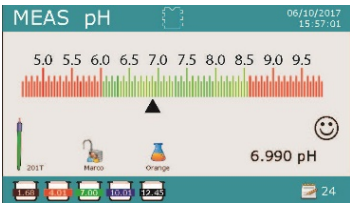
• **Simultaneous pH/Cond**

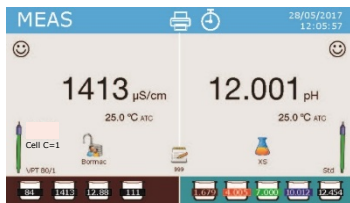Available only for PC8+DHS and PC80+DHS

- Connect the electrode to the BNC for pH/mV/ORP of device (yellow/blue).
- If the user does not use an electrode with incorporated temperature probe or an external NTC 30KΩ probe, it is advisable to manually update the temperature value (MTC).
- Remove the electrode from its protection cap, rinse with distilled water and gently dab with absorbent paper.
- Check the presence and eliminate any air bubbles in the bulb of the membrane by vertically shaking (as for the clinical thermometer). If it is present, open the side pad.
- Immerse the electrode in the sample while maintaining slight stirring.
- Consider truthful the measurement only when ☺ icon appears. To eliminate any errors due to user interpretation it is possible to use “HOLD” function (P6.8) which allows to freeze the measurement as soon as it achieves stability.
- **After the measurement, wash the electrode with distilled water and store it in the appropriate storage solution. Never store sensors in distilled water.**
- In the function “Display Complete P6.7”, the graphic representation of the electrode indicates the slope of the current calibration. 

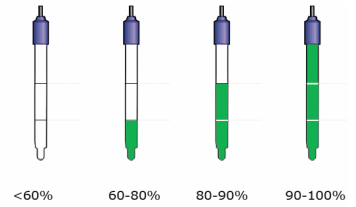• DHS Electrode

Electrodes equipped with DHS technology can save a calibration curve inside their memory. The calibrated sensor is automatically recognized by any DHS instrument and it acquires the calibration.

- Connect DHS electrode to the BNC and RCA of device.

- The device automatically recognizes the chip, the display shows information about the model and lot of the sensor and the date of the last calibration (if the electrode had already been calibrated).
- As soon as the DHS electrode is recognized the active calibration on the instrument becomes the one of the sensor (visible with the becker at the bottom left of the display or in P1.3 menu).
- If calibration is good (see the data report view and graphic representation -P1.3-) the electrode is ready to begin the measurements. Otherwise calibrate the electrode, the data will be automatically updated.
- A DHS electrode calibrated with an 8+DHS or 80+DHS device is ready to be used on any other pHmeter enabled DHS recognition and vice versa.
- When the electrode is disconnected, a message on the display informs the user of the deactivation of the sensor, the instrument retrieves its previous calibration, no data is lost!
- The DHS electrode does not require batteries and if used on a pH meter that is not enabled for chip recognition, it works like a normal analogue electrode.

- **Errors during calibration**

- **Error 1:** Not stable measure during calibration, before pressing  key wait for the stability .
- **Error 2:** The device does not recognize the buffers used for calibration.
- **Error 3:** The calibration has exceeded the limit time, only the points calibrated up to that time will be maintained.

10. ORP Parameter (Redox Potential)

pH8+DHS, pH80+DHS, PC8+DHS, PC80+DHS

- In measurement mode press  key to scroll through the parameters until you activate the screen **MEAS mV**.
- Connect the appropriate Redox to BNC for pH/mV/ORP (Yellow /Blue) measuring electrode and insert the sensor into the sample to be analyzed.
- Take the measure only when the stability icon  appears.
- ORP measurement is in mV and does not require a sensor calibration.
In order to verify the accuracy of the measurement, it is advisable to perform a quality control using a certified standard (200 / 475 or 650 mV).

ATTENTION: Before proceeding with the sensor calibration carefully consult the safety data sheets of the substances involved:

- Redox standard solutions.
- Storage solution for ORP electrodes.
- Filling solution for redox electrodes.

The careful reading of the safety data sheets of the solutions used, favors the elimination of residual risks related to skin contact, ingestion, inhalation or eye contact with the same which can generate, possible but not probable minor damage. Contact your local distributor for more information.

11. Measurement with Ion-Selective Electrodes (ISE/ION)

pH80+DHS, PC80+DHS

This series of devices can measure the concentration of ions as ammonium, fluorides, chlorides, nitrates etc. using a specific ion selective electrode for the ion of interest. Connect the electrode to the BNC pH / mV connector (yellow / light blue color).

- **Ion parameters setup**

- In measurement mode press “ MENU” key to enter into the SETUP menu.
- Use  and  keys to scroll until “Ion” menu and press  key to confirm.
- Use  and  keys to move and then select the program you want to edit.

The table below shows the structure of the setup menu for the Ion parameter, for each program there are the options that the user can choose and the default value.

Program	Description	Options	Factory default
P3.1	Measuring unit	mg/L – g/L – mol/L	mg/L
P3.2	Select low standard	0.001 ... 19999 ppm	0.001
P3.3	Set stability criteria	Stability / Time	Stability
P3.4	Last calibration Ion data	View / Print	View
P3.5	Due calibration Ion	NO – 1...99h – 1...99d	NO
P3.6	Ion Sensor ID	YES – NO	None
P3.7	Restore factory default	YES – NO	NO

P3.1 Measuring Unit

Access this menu to select the measuring unit for calibrating the instrument and making the sample reading.

- **mg/L** - factory default-
- **g/L**
- **mol/L**

Note: Use the same measuring unit in calibration and measurement.

If you change the unit of measurement, the calibration is automatically cancelled.

P3.2 Select Low Standard

Enter this menu to select the concentration of the first point of the calibration curve (more dilute standard). Other calibration points will be automatically identified by the software multiplying the concentration by factor of **10** (example: Low Standard 0.050 mg/L, the other calibration points will be 0.5 / 5 / 50 / 500 mg/L). The device can accept from a **minimum of 2** to a **maximum of 5** calibration points, when the second calibration point is terminated, the user can stop the calibration and save the points done up to that moment by pressing  key.

P3.3 Set Stability Criteria

Enter in this setup menu to set stability criteria to be used in calibration and in measurement mode.

- **Stability:** Equivalent to the "Medium" stability criterion for pH.
- **Seconds (0...180):** Use  and  keys to select after how many seconds to hold the measure (This is useful for volatile compounds). When using this option on the display, the countdown is activated, at the end of which the measurement is fixed.
- Press  to resume the time.

P3.4 Last Calibration Ion Data

Access this menu to get information about the last calibration performed.

- **“View”** (Factory default), a report is shown on the display with the following information about the current calibration:
CALIBRATION DATE and HOUR / TEMPERATURE / MEASURING UNIT / ID SENSOR / SLOPE %
- **“Print”**: connect the printer to the device through RS232 port (see outputs connections) to print the calibration report directly on paper format.

P3.5 Due calibration Ion

Enter this menu to set a calibration expiration, this option is essential in GLP protocol.

- No calibration expiration is set by default, use  and  keys to choose days or hours which must pass between two calibrations and confirm with  key.
- Once the calibration has expired, “” icon will appear on the display, anyway the user may still continue to make measurements.

P3.6 Ion Sensor ID

Enter this menu to assign an identifier name to the sensor currently in use.

- To activate this option, the external keyboard must be connected.
- The sensor ID will appear in the "Complete" display screen (P6.7) and in Complete or GLP paper print (P5.3)

P3.7 Restore factory default.

If the instrument does not work properly or incorrect calibrations have been performed press **Yes** with  key to return all Ion parameters to the factory default setting.

IMPORTANT: The factory reset of the parameters does not delete the stored data.

• Calibration with Ion-Selective Electrodes

Example for 2 points calibration 0.01 and 0.1 mg/L

- Enter in the **Ion** setup and select in program **3.1** the measuring unit **mg/L** and in the **3.2** the most diluted standard: **0.010**

The device automatically multiplies the lower standard for a factor by a 10 to locate the other points in the calibration curve. 


- Connect the appropriate ISE electrode to the connector for pH/mV/ORP (yellow/blue).

Important: If the ISE electrode is not combined, it is necessary to connect the specific reference electrode. Refer to the ISE electrode manual for any reference electrode filling electrolytes and for any ionic strength adjusters (ISA). 

- Press  key to get back to measurement mode and shift to **MEAS ION** with  key.
- Press and hold for 3 sec " CAL" key to enter into the calibration mode.

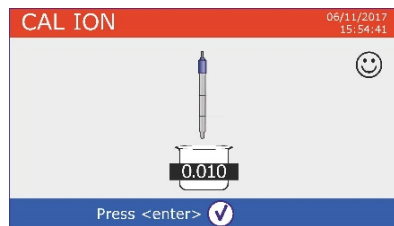- Rinse the electrode with distilled water  and gently dab with absorbent paper.

- Press  and immerse the electrode in the most diluted standard (Low Standard P3.2) as indicated by becker  0.010.


- When  icon appears (or when time runs out if you chose "Time" as the criterion for stability) press  key to confirm the first calibration point.


- Extract the electrode, rinse with distilled water  and gently dab with absorbent paper.

- Press  to continue and immerse the electrode in the next standard (**Low Standard X 10**) as indicated by becker  0.100.

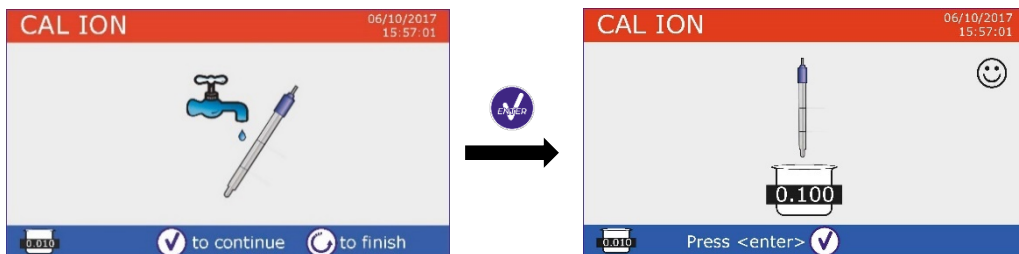- When  icon appears press  key to confirm the second calibration point.

- When the second point is finished, the user can exit calibration by pressing  key; instead press  to continue with the remaining points.
- When the calibration is completed, the calibration report is displayed with CALIBRATION DATE and HOUR / TEMPERATURE / MEASURING UNIT / ID SENSOR / SLOPE % for every range.

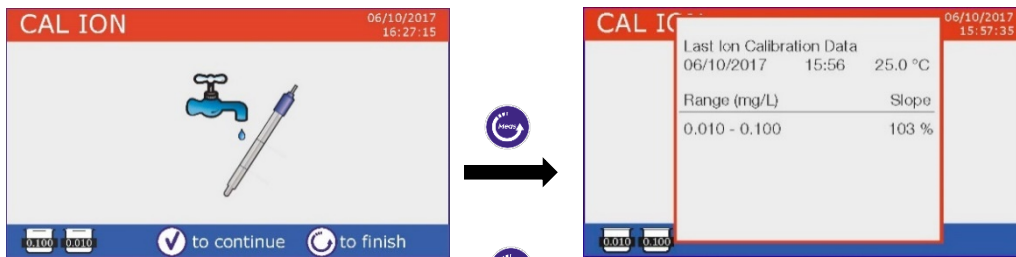Important: make at least two calibration points, if  key is pressed after the first calibration the display shows the error “Not enough calibrated points” and calibration is invalidated.


ATTENTION: Before proceeding with the calibration operations, carefully consult the safety data sheets of the substances involved:

- Standard calibration buffer solution.
- Storage solution for ISE electrodes.
- Filling solution for ISE electrodes.

• Measurement with Ion-Selective Electrodes

- Access the **Ion** setup menu to verify the correctness of the calibration and the instrumental parameters, Press  key to get back to in measurement mode and shift to key **MEAS ION** with  key.
- Correctly connect the ISE sensor to pH/mV/ORP BNC, rinse with distilled water, dab gently and insert it into the sample.
- The display shows the **measurement in mV** until the stability is reached.
- When the measure stabilizes, the mV measurement is replaced by the **concentration** of the analyte with the unit of measurement chosen by the user.

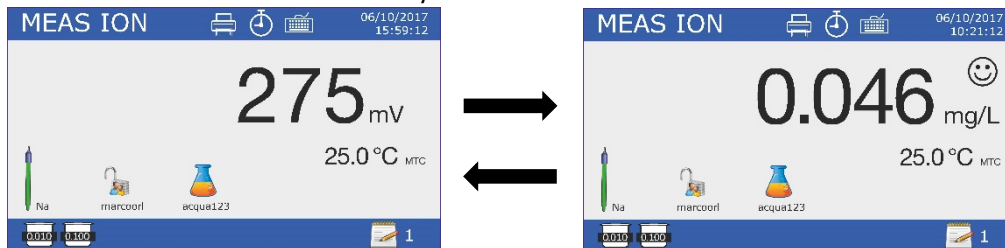Important: If the device is not calibrated, in measurement mode, only the mVs are displayed.

Note: If you are using as stability criteria countdown of seconds, to resume counting press .

12. Conductivity Parameter

Cond8+, Cond80+, PC8+DHS, PC80+DHS

Connect the Conductivity probe to the BNC type connector marked in grey, while the temperature probe must be connected to the RCA / CINCH Temp connector always on a grey background.

Conductivity is defined as the ability of the ions contained in a solution to conduct an electric current. This parameter provides a quick and reliable indication of the number of ions present in a solution.

• ... how to get to the conductivity?

The first Ohm law enunciates the direct proportionality in a conductor between the current intensity (I) and the applied electric potential difference (V) while the resistance (R) represents the constant of proportionality. In particular: “ $V = R \times I$ ”, consequently the resistance “ $R = V / I$ ”, where R = resistance (ohm), V = voltage (volt), I = current (ampere). The reverse of resistance is defined as conductance (G) “ $G = 1 / R$ ” and it is express in Siemens (S). The measurement of resistance or conductance requires a measuring cell, which consists in two poles of opposite charge. The reading depends on the geometry of the measuring cell, which is described through the cell constant parameter “ $C = d/A$ ” express in cm^{-1} where “d” is the distance in cm between the two electrodes and A their surface in cm^2 .

Conductance is transformed into a specific conductivity (k), which is independent of the cell configuration, multiplying it by the cell constant. $k = G \times C$ is expressed in S/cm even if measurement units “mS/cm” (1 S/cm \rightarrow 10³ mS/cm) and “ μ S/cm” (1 S/cm \rightarrow 10⁶ μ S/cm) are commonly used.

• Configuration of setup menu for Conductivity parameter

- In measurement mode press “ MENU” key to enter SETUP menu.
- Use  and  keys to scroll until “COND/TDS” menu and press  key to confirm.
- Use  and  keys to move and then select the program you want to edit.

The table below shows the structure of the setup menu for the COND/TDS parameter, for each program there are the options that the user can choose and the default value.

Program	Description	Options	Factory default
P2.2	Cell constant	0.1 - 1 - 10	1
P2.3	Calibration Solution	Standard / User	Standard
P2.4	Due calibration Cond	No – 1...99h – 1.99d	No
P2.5	Last calibration Cond data	View / Print	View
P2.6	Reference temperature for Cond	15...30 °C	25 °C
P2.7	Temperature compensation factor	0.0...10.0 %/°C – Ultrapure water*	1.91 %/C°
P2.8	Temperature calibration	YES – NO	No
P2.9	TDS factor	0.40...1.00	0.71
P2.10	Restore factory default	YES – NO	No

* Function available only for 80+ Series

P2.2 Cell Constant

Choosing the right conductivity cell is a decisive factor for obtaining accurate and reproducible measurements. One of the fundamental parameters to consider is to use a sensor with the correct cell constant in relation to the solution being analyzed.

The following table relates the sensor cell constant with the measuring range and the standard with which it is preferable to calibrate.

Cell Constant	0.1	1	10	
Conductivity Standard (25°)	84 μ S	1413 μ S	12.88 mS	111.8 mS
Measuring Range	0 – 200 μ S	200 – 2000 μ S	2 – 20 mS	20 – f.s. mS

Enter this setup menu to set the cell constant for the sensor you are using:

- **0.1**
- **1** -default-
- **10**

The cell constant in use appears on the lower left display. For each of the 3 selectable cell constants the device stores the calibrated points. Then, selecting the cell constant, the calibration points previously executed are automatically recalled

P2.3 Calibration Solution

Enter this setup menu to set the Automatic or Manual buffer conductivity standard calibration.

- **Standard:** -default- the device automatically recognizes max. 4 of the following standard **84 μ S/cm, 1413 μ S/cm, 12.88 mS/cm and 111.8 mS/cm.**
- **User:** the device can calibrate with one point defined by user.

Note: To obtain accurate results it is advisable to calibrate the device with standards close to the theoretical value of the solution to analyze.

P2.4 Due calibration Cond

Enter this menu to set a calibration expiration, this option is essential in GLP protocol.

- No calibration expiration is set by default, use  and  keys to choose days or hours which must pass between two calibrations and confirm with  key.

- Once the calibration has expired, “  ” icon will appear on the display, anyway the user may continue to make measurements.

P2.5 Last Calibration Cond Data

Access this menu to get information about the last calibration performed.

- “View” (Factory default): a report is shown on the display with the following information about the current calibration: CAL DATE and HOUR / TEMP / MEASURING UNIT / ID SENSOR / CELL CONSTANT for every range.
- “Print”: connect the printer to the device through RS232 port (see outputs connections) to print the calibration report directly on paper format.

P2.6 and P2.7 Temperature compensation in conductivity measurement it is not to be confused with the temperature compensation for pH measurement.

- In a conductivity measurement the value show on the display is the conductivity calculated at the reference temperature. So, it is corrected by the effect of temperature on the sample.
- Instead, in a pH measurement the value show is the pH at the displayed temperature. Here the temperature compensation involves the adaptation of the slope and the electrode offset at the measured temperature.

P2.6 Reference temperature for Cond

The conductivity measurement is strongly temperature dependent.

If the temperature of a sample increases, its viscosity decreases and this leads to an increase in ion mobility and of the measured conductivity, despite the concentration remains constant.

For each conductivity measurement the temperature to which it is related must be specified, otherwise it is a worthless result. Generally, the temperature refers to 25 °C or more rarely to 20 °C.

This device measures the conductivity at the actual temperature (ATC or MTC) and then convert it to the reference temperature using the correction algorithm chosen in P2.7 program.

- Enter this setup menu to set the temperature to which the conductivity measurement is to be reported.
- The device is able to refer the conductivity from 15 to 30 °C. As factory default setting is 25°C which is fine for most of the analysis.

P2.7 Temperature compensation factor

It is important to know the dependence on temperature (percentage variation of conductivity every °C) of the sample. To simplify the complex relationship between conductivity, temperature and ionic concentration different compensation methods can be used

- Linear Coefficient 0.00...10.0 %/°C** -default 1.91 %/°C - For the compensation of medium and high conductivity solutions, linear compensation can be used. The default factory value is fine for most of the routine measures.

Compensation coefficients for special solutions and for groups of substances are shown in the following table.

Sample	(%/°C)	Sample	(%/°C)
NaCl Saline Solution	2.12	1.5% Hydrofluoric acid	7.20
5% NaOH Solution	1.72	Acids	0.9 - 1.60
Diluted ammonia solution	1.88	Bases	1.7 – 2.2
10% Hydrochloric acid solution	1.32	Salts	2.2 - 3.0
5% Solforic acid solution	0.96	Drinking water	2.0

Compensation coefficients for calibration standards at different temperatures for $T_{ref} 25\text{ °C}$ are shown in the following table.

	0.001 mol/L KCl (147µS)	0.01 mol/L KCl (1413 µS)	0.1 mol/L KCl (12.88 mS)
0	1.81	1.81	1.78
15	1.92	1.91	1.88
35	2.04	2.02	2.03
45	2.08	2.06	2.02
100	2.27	2.22	2.14

The following formula is used to determine the calibration coefficient of a particular solution:

$$tc = 100x \frac{C_{T2} - C_{T1}}{C_{T1}(T_2 - 25) - C_{T2}(T_1 - 25)}$$

Where tc is the temperature coefficient to be calculated, C_{T1} and C_{T2} are conductivity at temperature 1 ($T1$) and temperature 2 ($T2$).

Each result with the compensated temperature is affected by an error caused by the temperature coefficient. The better the temperature correction, the lower the error. The only way to eliminate this error is to not use the correction factor, acting directly on the sample temperature.

Set as temperature coefficient 0.00% to deactivate the compensation. The displayed conductivity value refers to the real temperature value measured by the probe and not related to a reference temperature.

- **Ultrapure water***: Select this option when working with conductivity lower than 10 µS/cm. An icon in the upper left corner informs the user that this compensation mode is being used. When this threshold is exceeded, this option is automatically disabled, and linear compensation is activated.

* function available only for 80+ Series.

The temperature coefficient in ultrapure water varies strongly. The main reason for this is that the self-ionization of water molecules is more temperature-dependent than the conductivity caused by the other ions.

Note: Low conductivity measurements (<10 µS/cm) are strongly influenced by atmospheric carbon dioxide. To obtain reliable results it is important to prevent contact between the sample and the air, this can be achieved by using a flow cell or chemically inert gases such as nitrogen or helium that isolate the sample surface. 

P2.8 Temperature calibration

All the instruments within this series are pre-calibrated for a correct measure of the temperature. However, if there is a difference between the measured and the real one (usually due to a malfunction of the probe), it is possible to adjust the offset of + 5 ° C.

- Use  and  keys to correct the temperature offset value and press  key to confirm.

P2.9 TDS Factor

Access this menu to adjust TDS factor **0.4...1.00** / -default 0.71- to perform the conversion from conductivity to TDS.

- Consult paragraph “Other measurements perform with the conductivity cell”.

P2.10 Restore factory default.

If the instrument does not work properly or incorrect calibrations have been performed press **Yes** with  key to return all Cond parameters to the factory default setting.




IMPORTANT: the factory reset of the parameters does not delete stored data.

- **Automatic Conductivity calibration**

Example for one point calibration (1413 µS/cm) using a 1 cell constant sensor

- In **Cond** measurement mode press and hold for 3 sec “ CAL” key to enter in calibration mode.
- In PC devices Cond calibration can be accessed also from simultaneous pH/Cond screen then selecting **Cond**.

- Rinse the sensor with distilled water  and gently dab with absorbent paper. Flush with some mL of standard solution.
- Press  and immerse the sensor in 1413 $\mu\text{S}/\text{cm}$ standard solutions keeping it slightly in moving and making sure there are no air bubbles in the cell. In the becker all the Conductivity values flow that the device is able to recognize.
- When the value on the becker stabilizes on 1413 value and the  icon appears confirm by pressing  key.
- The display flashes the value actually measured, subsequently the calibration report is displayed showing the cell constant for each range and finally the device automatically returns to measurement mode. The icon becker relative to the calibrated point appears on the bottom left of display .
- *One point calibration is sufficient if the measurements are performed within the measurement range. For example, the standard solution 1413 $\mu\text{S}/\text{cm}$ is suitable for measurements between approx 200 - 2000 $\mu\text{S}/\text{cm}$.*
- To perform a multipoint calibration, once back in measuring repeat all calibration step. The becker icon related to the new calibrated point will be next to the previous one. It is advisable to start the calibration from the less concentrated standard solution and then continue in order of increasing concentration.
- When a new calibration of a previously calibrated point is carried out, it is overwritten on the old one and the cell constant is updated.
- For each cell constant (P2.2) the device memorises the calibration, in order to allow the user who uses multiple sensors with different constants to not be forced to recalibrate each time.
- The instrument recalls the last calibration relative to the parameters P2.2 (cell constant) and P2.3 (type of calibration solutions) selected.

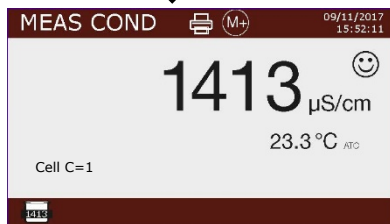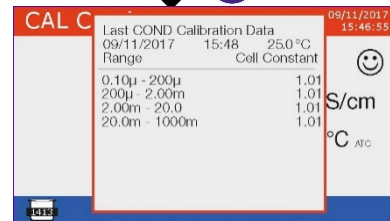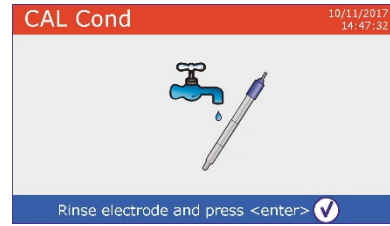Important: Standard conductivity solutions are more vulnerable to contamination, dilution and direct influence of CO₂ compared to pH buffers, which instead thanks to their buffer capacity tend to be more resistant. Furthermore, a slight temperature change, if not adequately compensated, can have significant effects on accuracy. Pay particular attention in the process of calibration of the conductivity cell and then be able to obtain accurate measurements. 

Important: To avoid contamination, always rinse the cell with distilled water before calibration and when switching from one standard solution to another to avoid contamination. Replace the standard solutions frequently, especially those with low conductivity. 

Contaminated or expired solutions may affect the accuracy and precision of the measurement.

ATTENTION: Before proceeding with the calibration operations, carefully consult the safety data sheets of the substances involved:

- Standard solution for calibration.


The careful reading of the safety data sheets of the solutions used favors the elimination of residual risks related to skin, contact, ingestion, inhalation or eye contact that can generate possible but not probable minor damages. Contact your local distributor for more information. 

• Manual conductivity calibration


Example for 5.00 $\mu\text{S}/\text{cm}$ calibration with 0.1 cell constant sensor

- Access **Cond/TDS** setup menu and set in P2.2 → **0.1** and P2.3 → **User**, press  key to return in measurement mode and switch on Cond mode.

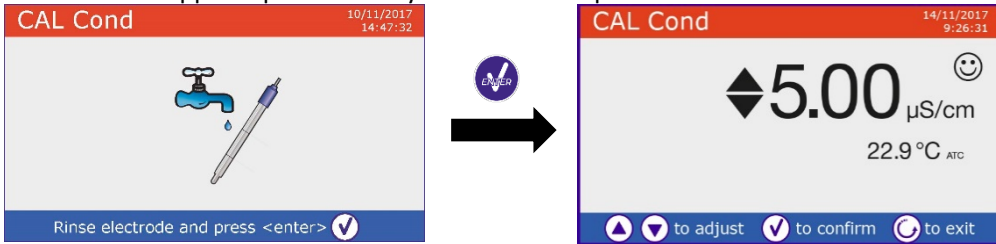- Press and hold for 3 sec “CAL” key to enter in calibration mode.
- Rinse the sensor with distilled water  and gently dab with absorbent paper.
Rinse with some mL of standard solution, press  and immerse the sensor in 5.00 µS/cm standard solution.
- Wait until the Conductivity value on the display stabilizes, when  icon appears use  and  keys to adjust the value inserting that of the standard solution (ex 5.00 µS/cm).
- When the  icon reappears press  key to confirm the point.

- Automatically the calibration report is displayed showing the cell constant for the whole range, press  to returns in measurement mode.
- The icon becker  for calibration is displayed at the bottom left, the value is preceded by the "U" letter indicating that the value **has been entered manually**.

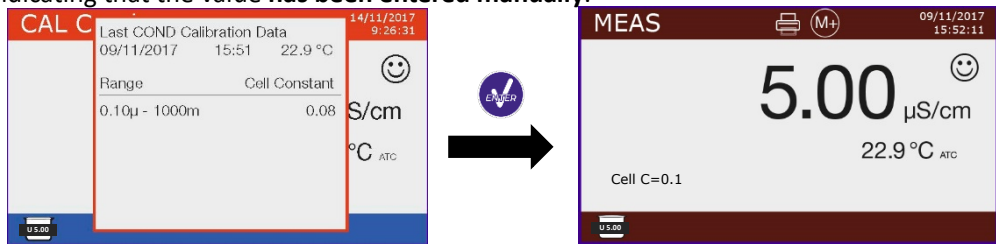- For each cell constant (P2.2) the device memorises the calibration, in order to allow the user who uses multiple sensors with different constants to not be forced to recalibrate each time.
- The instrument recalls the last calibration relative to the parameters P2.2 (cell constant) and P2.3 (type of calibration solutions) selected.

Note: if it is not aware of the exact coefficient of compensation, to get a calibration and an accurate measurement set in P2.7 → 0.00 % and then work bringing the solutions exactly to the reference temperature. Another method to work without temperature compensation is to use the appropriate thermal tables present on the majority of conductivity solutions.

Important: Always rinse the cell with distilled water before calibration and when switching from one standard solution to another to avoid contamination. Replace the standard solutions frequently, especially those with low conductivity. 


Contaminated or expired solutions can affect the accuracy and precision of the measurement.

• **Performing a Conductivity measurement**


- Enter the Conductivity menu setup to check the calibration and if necessary, update the reading parameters (see “Configuration of setup menu for Conductivity parameter Cond/TDS” section), press  key to return in measurement mode.
- Press  key to scroll through the different parameter screens until you find **MEAS Cond** (see paragraph “Operation of the device”).

The Cond measurement can be displayed on two different screens:

• **Digital**

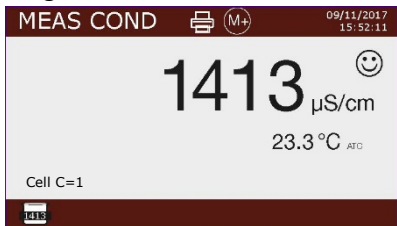• **Simultaneous pH/Cond**

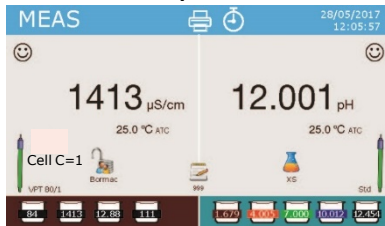Available only for PC8+DHS and PC80+DHS

- Connect the conductivity cell to the BNC for Cond of device (grey).
- If the user does not use a sensor with incorporated temperature probe or an external NTC 30KΩ probe, it is advisable to manually update the temperature value (MTC).
- Remove the electrode from its protection cap, rinse with distilled water and gently dab with absorbent paper **taking care to not scratch the electrodes**.
- Immerse the electrode in the sample, the measuring cell and any vent holes must be completely submerged.
- Keep slightly stirring, gently shaking the sensor to eliminate any air bubbles that would distort the measurement.
- Take the measurement only when 😊 icon appears. To eliminate any errors due to user interpretation, it is possible to use "HOLD" function (P6.8) which allows to freeze the measurement as soon as it achieves stability.
- The device uses six different measurement scales and two measurement units (µS/cm and mS/cm) depending on the value, the transition is carried out automatically by the device.
- Once completed the measurement, wash the conductivity cell with distilled water.
- The conductivity sensor does not require much maintenance, the main aspect is to make sure that the cell is clean. The sensor should be rinsed with plentiful distilled water after each analysis; if it has been used with water-insoluble samples previously clean it by immersing it in ethanol or acetone. Never clean the sensor mechanically, this would damage the electrodes compromising functionality. For short periods, store the cell in distilled water, while for long periods keep it dry.

13. Other measurements perform with conductivity cell

The conductivity measurement can be converted in TDS, Salinity and Resistivity parameters.

- In measurement mode press  key to scroll the different parameters TDS -> Salinity -> Resistivity.
- These parameters use conductivity calibration; pressing " CAL" key it will access directly to the calibration of conductivity.

• **TDS**

Cond8+, Cond80+, PC8+DHS, PC80+DHS

Total Dissolved Solids (TDS) correspond to the total weight of solids (cations, anions and non-dissociated substances) in one liter of water. Traditionally, the TDSs are determined by gravimetric method, but a simpler and faster method is to measure conductivity and then convert it to TDS by multiplying by the TDS conversion factor. Enter the P2.9 setup menu to modify the conductivity / TDS conversion factor.

Below are TDS factors in relation to the conductivity value.

Solution Conductivity	TDS Fcator
1-100 µS/cm	0.60
100 – 1000 µS/cm	0.71
1 – 10 mS/cm	0.81
10 – 200 mS/cm	0.94

The measurement of TDS is expressed in mg/L or g/L depending on the value.

• **Salinity**

Cond8+, Cond80+, PC8+DHS, PC80+DHS

Usually the UNESCO 1978 definition is used for this parameter, which provides the use of psu (Practical Salinity Units) as unit of measurement corresponding to the ratio between the conductivity of a sample of sea water and that of a standard KCl solution made by 32,4356 grams of salt dissolved in 1 kg of solution at 15 ° C. Relationships are dimensionless and 35 psu are equivalent to 35 grams of salt per kilogram of solution. Approximately 1 psu is equivalent to 1g / L and considering the density of water is equivalent to 1 ppt. The UNESCO 1966b definition can also be used which provides that the salinity in ppt is expressed with the following formula: $S_{ppt} = -0.08996 + 28.2929729R + 12.80832R^2 - 10.67869R^3 + 5.98624R^4 - 1.32311R^5$ Where R= Cond value (at 15°) / 42.914 mS/cm (Conductivity of Copenhagen Seawater Standard).

Note: If you want to perform a Salinity measurement in low conductivity solutions, it is advisable to deactivate the **Ultrapure Water** - P2.7 - temperature compensation.

• **Resistivity**

Cond80+, PC80+DHS

For low conductivity measurements such as ultrapure water or organic solvents it is preferable the resistivity value compared to the conductivity value. Resistivity is the reciprocal of conductivity $\rho = 1/\kappa$ (MΩ*cm).

14.Data Logger and printer functions

These series of devices are able to record values on internal memory or PC or to print them directly in different formats using an external printer to connect through the appropriate RS232 port. Recordings can be acquired manually or at preset frequencies.

To print: buy the printer separately, connect it with the RS232 cable to the connector labeled "Printer" on the rear panel of the instrument, connect the power supply to the mains and switch it on by pressing the switch **I / O**. For more information, see the technical manual of the printer. Use only the printer specified by the manufacturer. Contact your local distributor for purchase or further information. The original printer is already equipped with a roll of plain paper, power supply and RS232 cable for instrument connection.

PC connection: connect the USB type B cable inside each package to the USB port identified as "PC" on the rear panel of the instrument and the other end to a COM port on the computer.

Use only the USB cable supplied with each instrument.

• **Setup for Data Logger parameter**

- In measurement mode press “ MENU” key to enter the SETUP menu.
- Use and keys to scroll until “Data Logger” menu and press key to confirm.
- Use and keys to move and then select the program you want to edit.

The table below shows the structure of the setup menu for the Data Logger parameter, for each program there are the options that the user can choose and the default value.

Program	Description	Options	Factory default
P5.1	Save data	Memory – Printer – PC	Mem
P5.2	Data logging type	Man – Seconds – Minutes - Hours	Man
P5.3	Print format	Simple – Complete – GLP*	Simple
P5.4	Delete data in memory	YES - NO	No

* function available only for 80+ Series

P5.1 Save data

Enter this menu to select the location to save the recorded values:

- **Memory** - factory default - The recorded data is saved in the internal memory of the device. This serie of tools can store up to 1000 data in total, the progressive registration number appears alongside the icon .
- **Printer** Data is directly printed on an external printer connected through the RS232 port (see outputs connections). Enter P5.3 menu to select which information print in the header.

- **PC** The recorded data is directly sent to the PC through USB connection and processed using the appropriate software DataLink+.

P5.2 Data logging type

Enter this menu to select the acquisition mode data.

- **Manual** : The data is acquired or printed only when the user presses the key.
- **Seconds / Minutes / Hours** : By using the arrow keys set an automatic data acquisition frequency range. Press key to start and end the registration.

P5.3 Print format

Enter this menu to select which information to print in the header in a data logger printer mode.

- **For 8+ Series:**
 - **Simple:** progressive registration number / date and hour / value, UM and temperature.
 - **Complete:** device model / report last cal / progressive registration number / date and hour / value, UM and temperature.
- **For 80+ Series:**
 - **Simple:** progressive registration number / ID sample / date and hour / value, UM and temperature.
 - **Complete:** model and sn of the device / date last cal / IDs / progressive registration number / date and hour / value, UM and temperature.
 - **GLP:** model and sn of the device/ report last cal / IDs / progressive registration number / date and hour / value, UM and temperature.

Simple			Complete			GLP				
#1	16/02/2016	17:11:23	Model: PC80	16/02/2016	17:11:23	Model: PC80	16/02/2016	17:11:23		
Orange	5.502 pH	15.5 °C	Serial number: 162890220	Orange	5.502 pH	15.5 °C	Serial number: 162890220	Orange	5.502 pH	15.5 °C
	980 μS	15.5 °C	pH sensor ID: Standard	980 μS	15.5 °C	980 μS	15.5 °C	980 μS	15.5 °C	15.5 °C
#2	16/02/2016	17:11:33	Cond sensor ID: VPT801	#2	16/02/2016	17:11:33	Cond sensor ID: VPT801	#2	16/02/2016	17:11:33
Orange	5.512 pH	15.5 °C	Last Cal. pH: 2016-02-16 16:07	Orange	5.512 pH	15.5 °C	Last Cal. pH: 2016-02-16 16:07	Orange	5.512 pH	15.5 °C
	980 μS	15.5 °C	Last Cal. Cond: 2016-02-16 12:11	980 μS	15.5 °C	980 μS	15.5 °C	980 μS	15.5 °C	980 μS
			Company Name: XSinstruments.com							
			User ID: Administrator							
#3	16/02/2016	17:11:43		#3	16/02/2016	17:11:43		#3	16/02/2016	17:11:43
Orange	5.515 pH	15.5 °C		Orange	5.515 pH	15.5 °C		Orange	5.515 pH	15.5 °C
	980 μS	15.5 °C		980 μS	15.5 °C		980 μS	980 μS	15.5 °C	15.5 °C

• **Example 1 Data Logger mode**

Example for automatic pH recording on internal memory every 2 minutes.

- Enter in **Data Logger** setup mode.
- In **P5.1** select **Memory**.
- In **P5.2** select **Minutes**, with directional keys set “2” and confirm by pressing key.
- Press key to return in measurement mode and activate **MEAS pH** screen. In the upper band you can find icons - saving in internal memory - and - automatic data acquisition -.
- Press to start the registration, icon starts blinking, meaning that the memorisation is in progress. The number beside icon indicates how many data has been saved.
- Press again key to end the recording.

Note: automatic recording is suspended when you change the measurement parameter or enter the setup menu.

Note2: recording and printing in simultaneous pH/Cond screen occur for both the values.

• **Example 2 Data Logger mode**

Example for manual printing of Cond value with date of last calibration.

- Enter in **Data Logger** setup mode.
- In **P5.1** select **Printer**.
- In **P5.2** select **Man**.
- In **P5.3** select **Complete**.

- Press  key to return in measurement mode and activate **MEAS Cond** screen. In the upper band you can find icons  -print – and  - manual saving.
- Connect the printer to the RS232 port.
- When you want to print the value, press  key.

Note: For the first printed value of each series, the header selected in P5.3 is also printed.

• Recalling saved data on instrumental memory

- In measurement mode of the parameter of interest press  to display the data stored for that specific parameter.
- Using  and  keys scroll through the various pages of data.
- Press  key to exit and return in measurement mode.

Note: The first value of a series always has a progressive number "1" and is identified by an orange icon.

• How to delete data in memory

To delete data in instrumental memory, enter into P5.4 setup menu and select **Yes**.

Note: The stored data is kept even if accidentally the power supply fails.

IMPORTANT: The factory reset of the pH, Ion and Cond parameters does not delete the stored data.

15. Configuration menu

- In measurement mode press " MENU" key to enter into SETUP menu.
- Use  and  keys to scroll until "**Configuration**" menu and press  key to confirm.
- Use  and  keys to move and then select the program you want to edit.

The table below shows the structure of the setup menu for the **Configuration** parameter, for each program there are the options that the user can choose and the default value.

Program	Description	Options	Factory default
P6.1*	Sample ID	-	None
P6.2*	User ID	-	None
P6.3*	Company name	-	None
P6.4*	pH Sensor ID	-	None
P6.5*	Cond Sensor ID	-	None
P6.6	Password	-	None
P6.7*	Display information	Simple - Complete	Simple
P6.8	Reading with HOLD	YES – NO	NO
P6.9	Date format	dd/mm/yyyy – mm/dd/yyyy – yyyy/mm/dd	dd/mm/yyyy
P6.10	Date setting	-	-
P6.11	Time setting	-	-
P6.12	Temperature unit option	°C - °F	°C
P6.13	Set Language	Eng – Ita – Deu – Esp – Fra - Cze	English

* function available only for 80+ Series

P6.1 / P6.2 / P6.3 / P6.4 / P6.5 (function available only for 80+ Series)

To use these programs, connect the external keyboard to the USB port identified as "keyboard" (see paragraph "Input / output connection"). Use only the keyboard supplied by the manufacturer, for information contact the local distributor. The symbol  indicates that the connection has been successfully established. Access these programs to assign an alphanumeric identifier that will appear on the display in complete mode (P6.7) and during printing Complete and GLP (P5.3):

- **P6.1 Sample ID:** name of the sample; max 8 characters, appears on display along with the  icon and in printing for each measure in all formats.

- **P6.2 User ID:** name of the user; max 8 characters, appears in the print header Complete and GLP and on display along with the padlock: open  if no password is entered so the user can access all the menus, close  if password is entered so the user has access restrictions to the menus.
- **P6.3 Company name:** max 15 characters, appears only in the print header Complete and GLP.
- **P6.4 pH Sensor ID:** name of the pH sensor; max 8 characters, appears in the print header Complete and GLP and to display alongside the electrode calibration status .
- **P6.5 Cond Sensor ID:** name of the Cond sensor; max 8 characters, appears in the print header Complete and GLP and to display alongside the cell .

P6.6 Password

Enter this setup menu to insert, modify or disable the password.

- The **activated** password is indicated in the **Complete** display with the icon . If there are no limitations (Password **none**)  icon will appear.
- The user must provide the password **to calibrate the device, clear the saved data in data logger mode and and change date and time of the device.**

Note: In case of lost password, please contact Technical Assistance in order to unlock the instrument via Master Password which will be provided at the time.

For 8+ Series:

- 4-numeric characters password to enter with  and  keys, press  to move to the next number.
- To disable the password, enter as a new password “0 0 0 0”.

For 80+ Series:

- 4 alphanumeric characters password to enter exclusively with external keyboard.
- To disable the password, press “ENTER” as a new password with external keyboard.

P6.7 Display information *(function available only for 80+ Series)*

Enter this setup menu to select how many information to be shown on the display in measurement mode.

- **Simple** -default- The display only shows the measured value, the unit of measure, the temperature, stability criteria and cell constant.
- **Complete** On the display in measurement mode in addition to the default information, the various IDs and their icons are displayed (**P6.1 / P6.2 / P6.4 / P6.5**).

P6.8 Reading with HOLD

Enter this setup menu to enable or disable the HOLD stability criterion.

- **NO** -default- the measure is not fixed.
- **YES:** with this enabled option, the measurement is locked as soon as it reaches stability. The frozen value is indicated by the icon **HOLD**.

To unlock and resume the measurement until the next stability press .

P6.9 Date format

Enter this setup menu to modify the date format.

- **dd/mm/yyyy** -default-
- **mm/dd/yyyy**
- **yyyy/mm/dd**

P6.10 / P6.11 Date and Time Setting

Enter this setup menu to update device date and time.

P6.12 Temperature unit option

Enter this setup menu to select which temperature measurement unit to use.

- **°C** -default-
- **°F**

P6.13 Set Language

Enter this setup menu to select which language to use.

- English -default-
- Italiano
- Deutsch
- Espanol
- Francais
- Czech

16. Magnetic stirrer with independent control

All the device of 8+ and 80+ Series can be connected to the appropriate magnetic stirrer independently controlled; it is included if the instrument is purchased in the STIRRER version or it can be ordered separately. The magnetic stirrer is equipped with a functional probe holder with 3 compartments to hold the measurement sensors and a hole for the temperature probe. Use only the original stirrer. Contact your local distributor for more information.

• Functioning


- Connect the device (see “Outputs / Inputs connections”) to the socket on the back panel of the stirrer exclusively through the cable supplied.
- Press key to switch on the stirrer, the first green LED lights up.
- To activate the engine press once , then use and keys to adjust the speed. The stirring speed is adjustable on 15 levels, every 5 steps a green LED light up.

Note: It is possible to fix the stirrer both on the right and on the left side of the instrument by means of special metal foils on the base of the stirrer (example in the picture below). Unscrew the foils and then screw them with a hole under the device and the other under the stirrer.

The positioning made in this way eliminates the residual risk of possible slight damages caused by manual handling of loads or by contact with moving parts.

Controller of magnetic stirrer

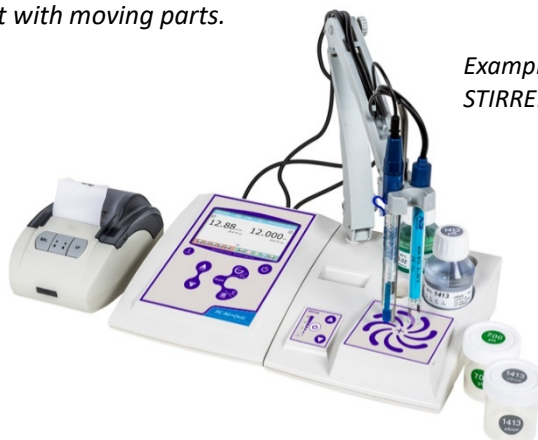Example of multiparameter STIRRER version

Compliance with the above, and reading of this manual before each use, eliminates the residual risk of minor damage from electrocution as far as, they are, however, unlikely.

17. Software DataLink+ (for Windows 7/8/XP/10)

It is possible to connect the Series 8+ and Series 80+ devices to PC and then use DataLink+ software to perform data download, datalogger directly to PC and Excel and PDF file exports with headers and ID.

You can download the software free on the websites (pay attention to the drivers to be installed):

- https://www.giorgiobormac.com/it/download-software_Download.htm
- https://www.xsinstruments.com/en/download_000034.htm
- Connect the device USB type B port with the appropriate USB cable supplied (see “Input/Output connection”) to a COM port of PC;
- Connect the USB type B cable inside each package to the USB port identified "PC" on the rear panel of the instrument and the other end to a COM port on the computer.

- Use only the USB cable supplied with each instrument.
- Start the program and then turn on the instrument.
- Wait for the connection to be established (the connection data are displayed at the bottom left of the display).

• Functions

- **Download:** the data stored in the instrument memory are downloaded to a PC and displayed in the table for processing;
- **M+:** instant acquisition of a value (is equivalent to manual datalogger);
- **Logger:** automatic acquisition with set frequency;
- **Clear:** emptying of the data in the table. If the password is active, it will be requested;
- **Export to Excel / Export to PDF:** export to PDF or to Excel of all the data in the table, of graphs, of calibration report and device information;
- **Save to file / Open from file:** saving the data in the table and possibility to reload them in order to process or continue recording;
- **Select language:** set the interface language (Eng – Ita – Deu – Esp – Fra – Cze);
- **Table / Chart:** display mode of the acquired data. The graphs are divided by parameter and can be printed separately.

The screenshot displays the software interface with several key components:

- Functions:** A menu bar at the top left contains icons for Download, Export to Excel, Export to PDF, Clear, Save to File, and Chart.
- Display of data acquired in the table or chart:** A large table in the center shows columns for Sample ID, #, Date, Time, Value, MU, Temp, MU, MIC/ATC, and DHS. The table contains 19 rows of data.
- Reproduction of the instrumental display for real-time measurements:** A small window in the top right corner shows real-time measurements for two parameters: 4,18 and 4,163.
- Calibration report for each parameters:** A window on the right side displays detailed calibration information, including constants, ranges, and buffers for various parameters.
- Automatic data-logger:** A status bar at the bottom left shows connection details: Connection established, Port COM4, Mode PCSD, Baud 9600, Parity N, Stop 8, Flow Control None, and a data logger icon.
- Connection data and device information:** A status bar at the bottom right shows system information: CPU 40, RAM 4448/1024, and a date/time display of 11:29.

18. Warranty

- The manufacturer of this device and its accessories offers the final consumer of the new device the three-year warranty from the date of purchase, in the event of state-of-the-art maintenance and use.
- During the warranty period, the manufacturer will repair or replace defective components.
- This warranty is valid only and exclusively on the electronic parts of the device and does not apply, if the product has been damaged, used incorrectly, exposed to radiation or corrosive substances, if foreign materials have penetrated inside the product or if changes have been made, which have not been authorized by the manufacturer.

19. Disposal of electronic devices

This equipment is subject to the regulations for electronic devices.
Dispose of in accordance with local regulations.

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