

Centrifuges made by Funke-Gerber
approved - exact - reliable

**FUNKE
GERBER**

SuperVario N

Multi-Purpose Centrifuge for the Dairy Industry



Multi-Purpose Centrifuge SuperVario N

Gerber-Centrifuge up to 36 Butyrometers, Babcock fat determination, Röse-Gottlieb fat determination* and solubility index determination (ATPI)

*Pay attention to existing safety rules!

Features:

- **Easy Handling:**
5 selectable modes for centrifugation
 (Gerber test, Röse-Gottlieb, Babcock, Solubility and Free Mode)
Free Set-up (Free Mode) of
 1. **Speed/RCA**** (adjustable from 600 to 1200 rpm, the corresponding RCA is shown simultaneously)
 **RCA = Relative Centrifugal Acceleration
 2. **Temperature** (adjustable from amb. temp. to 68°C)
 3. **Centrifugation time** (adjustable from 1 to 99 minutes)
- **Centrifugation drum of stainless steel**
- **Security lid locking device**
- **Automatic unbalance switch off**
- **Electrical brake**

Technical Specification:

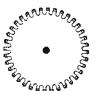
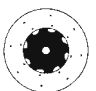

Connection: 230V/50...60 Hz/1200 VA
 Weight (empty): 26 kg
 Full height: 460 mm
 Operation height: 370 mm
 Speed range: 600...1200 UpM
 Temperature range: Amb. temperature to 68°C



Order Data:

Name	Art. No.
Centrifuge SuperVario-N (without heads)	3680

Accessories:

Name	Art. No.	
 A	Head A (up to 36 Butyrometer/up to 18 Babcock-bottles)	3685
	Butyrometer-bucket of cast aluminium	3631
	Babcock-bucket	3632
 B	Head B (Protective drum for 8 Mojonnier tubes)	3686
 C	Head C (up to 6 solubility index tubes)	3687
	Solubility index tube bucket	3633

Distribué par :

Z.A de Gesvrine - 4 rue Képler - B.P.4125
 44241 La Chapelle-sur-Erdre Cedex - France
 t. : +33 (0)2 40 93 53 53 | f. : +33 (0)2 40 93 41 00
 commercial@humeau.com

