



## CalCheck | LED Calibration Verification


CalCheck is a reusable LED device for quick and reliable luminometer calibration verification – in less than one minute!

Though Hygiena™ luminometers already run a calibration self-check at start up, a robust quality control program often requires proof of instrument calibration. In-house calibration checks confirm the Hygiena luminometer is working properly and demonstrates due diligence that a monitoring system is in control. CalCheck provides an all-in-one, reusable positive and negative calibration verification device, activated by the click of a button.

### Benefits:

- Verifies instrument calibration in-house as part of a quality control program
- Demonstrates due diligence to auditors and proves that your system is working properly
- One time purchase allows for unlimited calibration verification
- Powered by an easily replaceable coin battery
- Stable LED-based light source
- Compatible with Hygiena EnSURE™ and SystemSURE Plus luminometers

### CalCheck Procedure

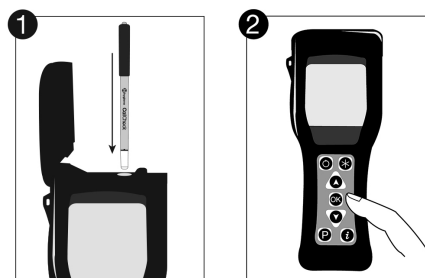

Catalog No.

Description

Quantity

CAL

CalCheck LED Calibration Verification Device

1



Find support documents, instructional videos, and more at [www.hygiena.com](http://www.hygiena.com)